



Making and Restoring Backups using Wi-Fi

Proloquo2Go (iPad, iPhone and iPod touch). In this tutorial you will learn how to make and restore backups with your web browser using Wi-Fi.

Turning on Wi-Fi



Tap on the Home button on your device and navigate to the Settings App (it looks like a gear). Look for the Wi-Fi section and tap on Wi-Fi.



Make sure WiFi is set to ON **1** and that you connect to the same Wi-Fi network on which your computer is located. If you are not currently connected to that network, tap the name of the network and enter the password if necessary. Then tap on the Home button on your device.



Launch Proloquo2Go by tapping on its icon. Go to the Options view.

*On an iPod touch or iPhone, the Options View can be accessed by tapping the views toolbar button at the bottom left of the screen, and then tapping the Options button **1** in the Views display.*



*On iPad, the Options button **1** is the right-most button on the toolbar.*



Tap Backup and then tap Backup Now. **1** This will make a full backup of your current settings and vocabulary customizations.


Transferring a Backup to the Computer




Tap on Connect to Computer. **1** Next, read the instructions on the subsequent screen and press Start. Then read the instructions on the next screen.




Do not exit this screen until after you are done transferring backups. While downloading backups the device should not go to sleep mode as this disrupts the connection between the app and the computer. It helps to keep your device connected to a power source.

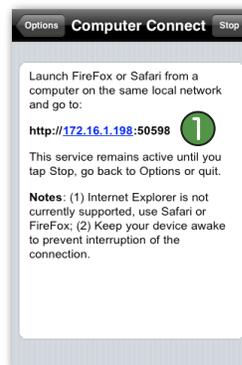


Launch Safari or FireFox web browser on your computer. (note that Internet Explorer is not currently supported), and enter the URL that is displayed in the Proloquo2Go Computer Connect screen. **1**

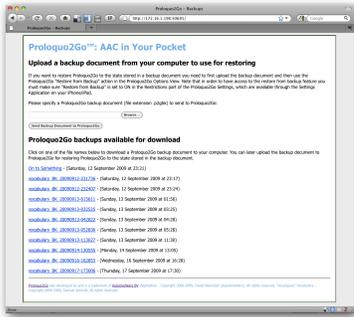

Press enter on the computer after typing the URL so that your browser will try to connect to Proloquo2Go. You will see a screen similar to the one on the left.

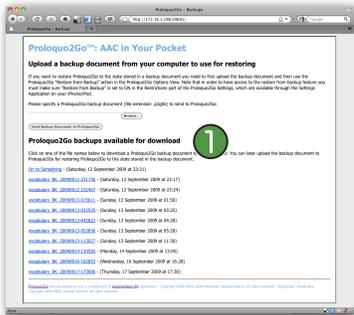

To download a backup to your computer, simply click on one of the backups listed under "Proloquo2Go backups available for download".

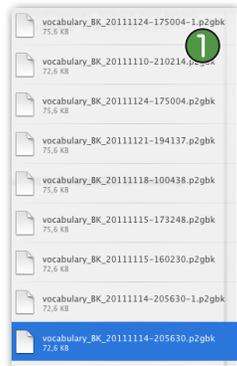


The backup will be downloaded as a .p2gbk file. Note that while you may change the name of the download backup file, you should not remove or change the file extension (it should remain .p2gbk).


Transferring a backup from your computer to Proloquo2Go




If you are not in the Computer Connect screen in Proloquo2Go anymore, go back to the Connect to Computer screen in the Proloquo2Go Options and then follow steps 1 through 5 of *Transferring a backup to the computer* above.

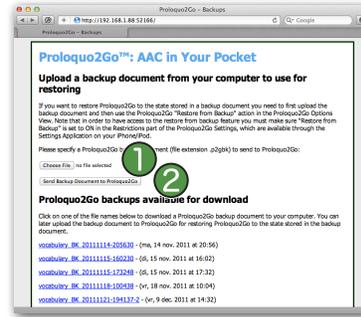

Click the Choose File button and select the backup file you want to send to Proloquo2Go. Now press the Send Backup document to Proloquo2Go button.

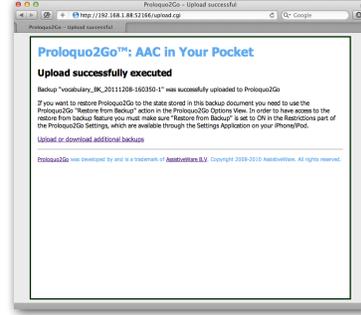

When the transfer is completed you will get a confirmation screen on your computer that the transfer was successful. Close the browser window in your computer.

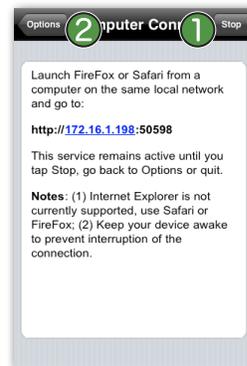

Tap the Stop button of the screen in Proloquo2Go. Then tap the Options button to go back to the Options.

Restore Proloquo2Go to a previous backup



In the Options screen, tap Backup and then tap "Restore From Backup" (if this item is not shown, tap Options on top of the Backup screen, tap Restrictions and toggle Restore from Backup to ON). Now choose the backup to restore from and tap the Restore from Backup button to confirm.